



SEPARATION

SINGLE-STAGE PROTEIN PURIFICATION

ROTOR STATOR HIGH-GRADIENT
MAGNETIC SEPARATOR MES-RS

ANDRITZ

ENGINEERED SUCCESS

Engineered for maximum yield

Biopharmaceutical active ingredients are proteins that have to be isolated from the fermentation broth in which they are produced. These purification processes always involve several process steps and each of these time-consuming and cost-intensive steps leads to losses of the target protein. Overall, this results in huge losses of the valuable component. The ANDRITZ rotor stator high-gradient magnetic separator MES-RS was developed to overcome these losses and extract more from less.

SINGLE-STAGE PROTEIN PURIFICATION

The “all-in-one” approach of the ANDRITZ MES-RS simplifies the conventional cascade of process steps, resulting in a single unit for the entire process. An impressive example of the efficiency of this technology is the purification process for equine chorionic gonadotropin (eCG). The standard purification process achieves an overall yield of 28%. When using the MES-RS, the protein harvesting step results in a yield of 97%, which means an overall process yield of 75%. Starting from the same batch size of fermentation broth, this means the final product output is three times higher. The magnetic separator enables the extraction of a single protein fraction directly from the non-purified complex feedstock, thus setting new standards in efficiency and yield intensification.

The MES-RS concept for protein purification, which utilizes high-gradient magnetic separation (HGMS) technology, allows ANDRITZ Separation to offer a novel downstream process for biopharmaceuticals. The state-of-the-art technology for protein purification is liquid chromatography, which requires multi-stage pretreatment to be able to process a complex feedstock. The MES-RS therefore reduces the effort for downstream processing while substantially increasing the yield.



Lab-scale MES-RS in operation, © KIT, Karlsruhe

MAIN APPLICATIONS


- Monoclonal antibodies (mAb)
- Hormones
- Enzymes
- Vaccines
- Biopharmaceutical proteins
- Feedstocks with extremely low titer
- Feedstocks with high downstream effort
- Magnetic particles

YOUR BENEFITS

- Avoids multi-stage pretreatment steps in downstream processing
- Single unit for prepurification, including first liquid chromatography step
- Target protein yield of up to 95%
- Purity of up to 97% after single-stage process
- Reduction of processing time to just several minutes
- Purification of broths that previously could not be purified
- Processing in accordance with GMP and ATEX requirements
- Fully automatic operation

PROCESSING PARAMETERS

Separator type	MES-RS 25	MES-RS 100	MES-RS 1000
Process chamber volume [l]	0.25	1	10
Loaded particle capacity [g]	25–100	50–500	500–5,000
Flow rate [l/min]	0.5	3	30
Fermenter batch size [l]	0.5–10	2–500	50–5,000
Operating principle	Highly selective binding of proteins to functionalized bead surface		
Separation principle	Magnetic capture of magnetized, functionalized particles		
Operating temperature	5–40 °C		
Operating mode	Batch		
Operating hours	24/7, fully automatic		
Construction material	Stainless steel		
Optional feature	Adjustable magnetic flux density		

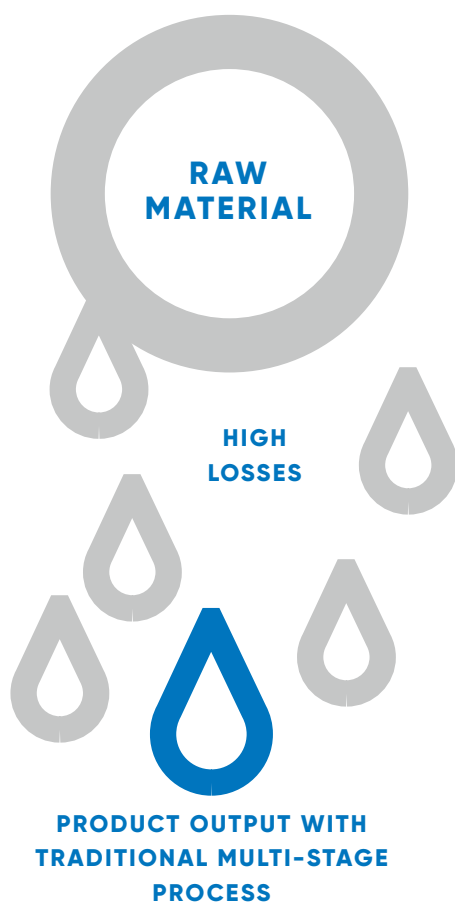

Traditional purification process: Very high losses and thereby poor efficiency in the multi-stage process



The ANDRITZ MES-RS process: Minimized losses in single-stage isolation process

Highly efficient extraction and purification of proteins through magnetic separation

OPERATING PRINCIPLE

In contrast to filtration or sedimentation techniques, the MES-RS separates only magnetizable components from the fermentation broth, and neither density nor particle size have any influence on the separation efficiency. The principle of magnetic separation, which is described in the following four steps, provides a highly selective method of extracting and purifying proteins.

Magnetic beads coated with a highly selective, functionalized surface (magnetic particles) are used to adsorb a specific protein fraction (step 1).

The process chamber of the MES-RS, which separates the magnetic particles, consists of an arrangement of multiple separation discs. The magnetic field magnetizes the particles and the separation discs, which causes the two components to attract (step 2).

By applying a magnetic field to the magnetic particles, it is possible to extract both the magnetic particles and the adsorbed proteins.

Since an electromagnet is used, the separation forces can be switched on and off. In addition, it is possible to rotate every second separation disc (rotor stator principle), which allows the captured magnetic particles to be reconverted into a slurry (step 3).

Both features of the MES – the electromagnet and the rotor stator design – make it possible to carry out several process steps within one process chamber. The concentrated and purified target protein can be discharged directly out of the MES-RS process chamber (step 4).

STEP 1:


BINDING OF TARGET PROTEIN

- 1 Mixing of magnetic particles and fermentation broth, binding of target protein to highly functionalized magnetic particle surface

STEP 2:

SEPARATION OF MAGNETIC PARTICLES


- 2 Feeding of magnetic slurry into the process chamber
- 3 Separation of particles on separation discs by magnetization
- 4 Discharging of particle-free waste slurry until the feed tank is empty

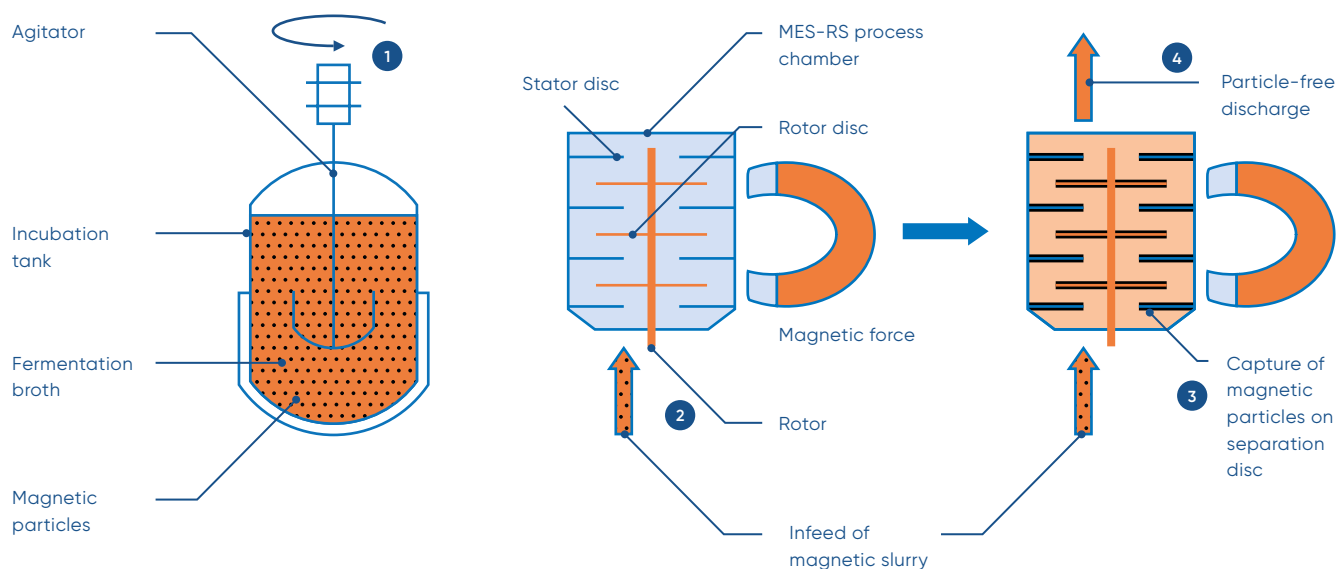

STEP 3:

WASHING OF PARTICLES


- 5 Capture of particles and displacement of mother liquor by wash liquor
- 6 Reconversion of particle cake into slurry to remove impurities from the particle cake
- 5 6 Repetition as often as necessary

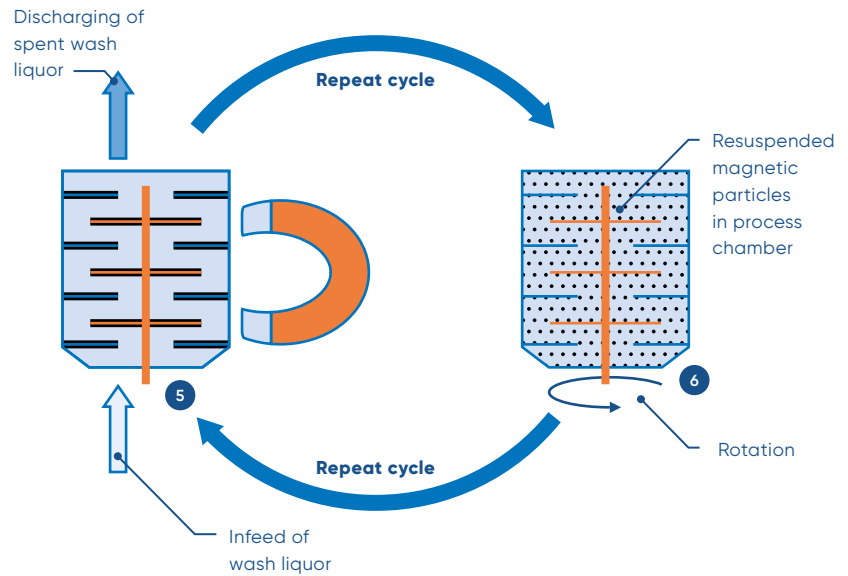

STEP 4:

ELUTION OF TARGET PROTEIN

- 7 Capture of particles and displacement of wash liquor by elution liquor until the process chamber is filled with 100% elution liquor
- 8 Reconversion of particle cake into slurry to elute target proteins from magnetic particles
- 9 Capture of particles and displacement of protein-enriched elution liquor
- 8 9 Repetition as often as necessary

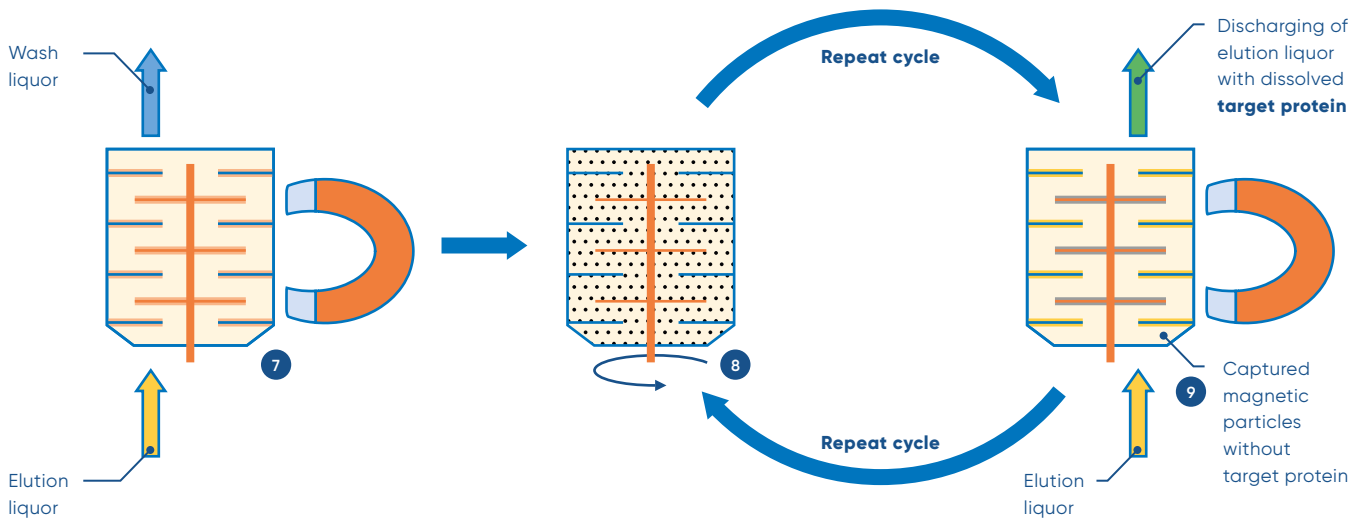

"The ANDRITZ magnetic separator is the next step in chromatographic process intensification and therefore has the potential to revolutionize product purification processes in biotechnology."

PROFESSOR DR. ING. MATTHIAS FRANZREB


Inventor and expert in biopurification

Technical data

Separator type	Process chamber volume [l]	Magnetic field strength [T]	Depth [mm]	Width [mm]	Height [mm]	Weight [kg]	Drive [kW]	Power [kW]
MES-RS 25	0.25	0.2	750	1,400	1,800	750	1	3.6
MES-RS 100	1	0.2	750	1,400	1,800	750	1	3.6
MES-RS 1000	10	0.2	1,000	1,800	1,800	1,500	3	15

All technical data are approximate and subject to change without notice.

Piping and instrumentation

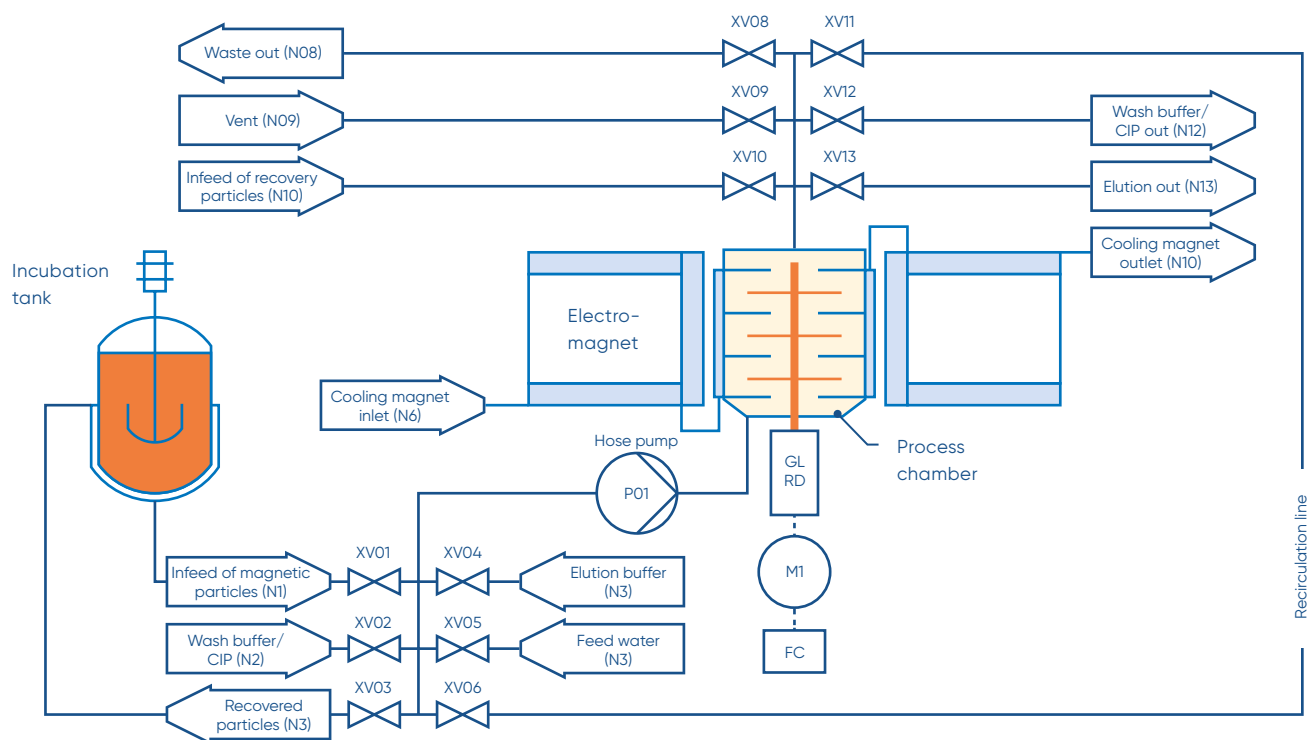

The piping and instrumentation diagram (P&ID) shows the setup of the MES-RS process skid. Depending on the process, the built-in actuators are used to carry out the process steps necessary for the purification process.

Your full-service provider

With ANDRITZ, you gain access to one of the world's largest OEM manufacturers for solid/liquid separation systems, including such well-known brands as 3Sys Technologies, Bird, Delkor Capital Equipment (Pty) Ltd., Escher Wyss dryers, Frautech, Guinard Centrifugation, KHD Humboldt Wedag, Krauss-Maffei centrifuges, dryers, and filters, Lenser, Netzsch Filtration, Rittershaus & Blecher, Royal GMF Gouda, Sprout Bauer, and Vandenbroek. Whether you need spare parts, rentals, local service, repairs, upgrades, or modernization of your equipment,

ANDRITZ is your true full-service provider. From initial consulting through to service agreements, process optimization, and training programs, we are always looking for ways to minimize downtime and increase predictability in operations while raising your overall production efficiency. Wherever you operate, our network of 550 service specialists and global service centers ensures we'll always be there to support you for many life cycles to come. Let's sit down and see how we could take your operations to the next level.



LOCAL SUPPORT



REPAIRS & UPGRADES



SECOND-HAND & RENTALS



TRAINING



OEM SPARE PARTS



SERVICE AGREEMENTS



PROCESS OPTIMIZATION



LAB AND ON-SITE TESTS

Intelligence for machine and process control

METRIS addIQ **CONTROL SYSTEMS**

With Metris addIQ, you have a well-proven, intelligent control solution for industrial processes and machines. Our solid/liquid separation specialists use their in-depth expertise to provide scalable solutions that are individually tailored to regional and application requirements. Whether you're automating new equipment or upgrading to extend the lifecycle of existing systems, we find the ideal solution for you.

Metris
ANDRITZ Digital Solutions

ANDRITZ offers a broad and constantly growing range of innovative products and services in the industrial digitalization sector under the brand name Metris, helping customers to enhance plant efficiency and profitability, optimize the use of resources, achieve constant and highest product quality, reduce production downtime, and maximize user-friendliness.



WHAT'S YOUR SEPARATION CHALLENGE?

ANDRITZ Separation is the world's leading separation specialist with the broadest technology portfolio and more than 2,000 specialists in 40 countries. For more than 150 years, we have been a driving force in the evolution of separation solutions and services for industries ranging from environment to food, chemicals, and mining & minerals. As the OEM for many of the world's leading brands, we have the solutions and services to transform your business to meet tomorrow's changing demands – wherever you are and whatever your separation challenge. **Ask your separation specialist!**

AFRICA

ANDRITZ Delkor (Pty) Ltd.
p: +27 11 012 7300
separation.za@andritz.com

EUROPE

ANDRITZ Separation GmbH
t: +49 8139 80299 0
separation.de@andritz.com

ASIA

ANDRITZ Singapore Pte. Ltd.
p: +65 6512 1800
separation.sg@andritz.com

NORTH AMERICA

ANDRITZ Separation Inc.
p: +1 817 465 5611
separation.us@andritz.com

AUSTRALIA

ANDRITZ Pty. Ltd.
p: +61 3 8773 4888
separation.au@andritz.com

SOUTH AMERICA

ANDRITZ Separation Ltda.
p: +55 47 3387 9100
separation.bra@andritz.com

CHINA

ANDRITZ (China) Ltd.
p: +86 757 8258 6802
separation.cn@andritz.com

ANDRITZ.COM/SEPARATION

JOIN US ON SOCIAL MEDIA



All data, information, statements, photographs, and graphic illustrations in this leaflet are without any obligation and raise no liabilities to or form part of any sales contracts of ANDRITZ AG or any affiliates for equipment and/or systems referred to herein. © ANDRITZ AG 2021. All rights reserved. No part of this copyrighted work may be reproduced, modified, or distributed in any form or by any means, or stored in any database or retrieval system, without the prior written permission of ANDRITZ AG or its affiliates. Any such unauthorized use for any purpose is a violation of the relevant copyright laws. ANDRITZ AG, Stattegger Strasse 18, 8045 Graz, Austria. PB Magnetic Separator MES-RS 2.0/04.2021 EN

